

Southern Illinois University Edwardsville

SPARK

SIUE Faculty Research, Scholarship, and Creative Activity

Spring 4-16-2020

Extensible Interface for a Compact Spectrophotometer for Teaching Molecular Absorption in the Undergraduate Laboratory

Edward Navarre
enavarr@siue.edu

Follow this and additional works at: https://spark.siu.edu/siue_fac

 Part of the [Analytical Chemistry Commons](#)

Recommended Citation

Navarre, Edward, "Extensible Interface for a Compact Spectrophotometer for Teaching Molecular Absorption in the Undergraduate Laboratory" (2020). *SIUE Faculty Research, Scholarship, and Creative Activity*. 141.

https://spark.siu.edu/siue_fac/141

This Article is brought to you for free and open access by SPARK. It has been accepted for inclusion in SIUE Faculty Research, Scholarship, and Creative Activity by an authorized administrator of SPARK. For more information, please contact magrase@siue.edu, tdvorak@siue.edu.

An extensible interface for a compact spectrophotometer for teaching molecular absorption in the undergraduate laboratory

Edward C. Navarre *

5 Department of Chemistry, Southern Illinois University Edwardsville, Edwardsville IL 62026-1652

ABSTRACT

13 A simple computer interface for controlling a compact spectrograph for use as a spectrophotometer in
14 an undergraduate teaching laboratory was developed. The project was implemented on a Raspberry Pi
15 computer which permits the integration of a light source into the software. The interface was written
16 in Python to facilitate modification by the user and because of its compatibility with several computing
17 platforms. An implementation of the project in Linux on a Raspberry Pi computer is described.
18
19
20
21
22
23
24

KEYWORDS

25 First-Year Undergraduate / General; Second-Year Undergraduate; Laboratory Instruction; Hands-On
26 Learning / Manipulatives; Spectroscopy; Laboratory Equipment / Apparatus
27
28
29

INTRODUCTION

30 15 Absorbance spectrophotometry is widely used throughout the chemistry curriculum to the
31 extent that it is described at least qualitatively in most introductory chemistry textbooks and often
32 basic univariate quantitative methods are included.^{1,2} The theoretical and quantitative aspects of
33 absorbance spectrophotometry are often presented for the first time in a classical quantitative
34 chemical analysis course^{3,4} taken in an undergraduate student's second or third year at university.⁵
35 As a result, simple and cost-effective equipment is essential to teaching this technique in a variety of
36 university laboratory situations with levels of student sophistication from first-semester university
37 students to students in their final semester of an undergraduate degree. Many commercial solutions
38 are offered ranging from about \$400 to \$2000 for basic absorbance systems and considerably more for
39 advanced features. Additionally, many publications in this journal describe low-cost approaches to
40 making absorbance spectrophotometry more accessible^{6,7} and innovative ways to use low-cost
41 instruments.^{8,9}
42
43
44
45
46
47
48
49
50
51
52
53
54
55
56
57
58

1
2
3 The author was recently faced with the decision of how to replace a set of Bausch and Lomb
4 Spectronic 20 used in chemistry teaching laboratories. Although the ThermoScientific Spectronic 200
5 would be a 'one-for-one' replacement, the author pursued a different premise: how could existing
6 spectrophotometric resources be used to solve the replacement issue? The author's academic
7 department possessed several compact spectrographs that were originally acquired for an upper-level
8 teaching laboratory. The technical solution described below is the author's attempt to develop a
9 compact, economical replacement to the functionality of the Spectronic 20 in a sophomore-level
10 quantitative analysis laboratory. The project sought to replicate the features without adding
11 complexity to the use of the instrument.
12
13
14
15
16
17
18
19
20

21 An additional motivation for the project was two limitations related to the manufacturer-
22 provided software: space considerations and complexity. The laboratory into which this project would
23 be deployed has limited space and an earlier trial of the compact spectrographs with small form factor
24 PCs proved too unwieldy for regular use. Additionally, the compact spectrograph manufacturer
25 justifiably wants to offer all of the hardware features in the software interface. This wealth of options
26 was more confusing than helpful to students in the early trial of the hardware.
27
28
29
30
31
32

33 The author pursued an interface design that presented all of the important spectrophotometric
34 conditions to the student in one place. As a result, the student user would be able to focus on
35 learning about spectrophotometry rather than spend time understanding a software interface.
36 Presenting only the essential experimental conditions limits the available spectrometer functions by
37 removing features like derivatives, smoothing, and event triggers. Additionally, convenience features,
38 such as automatic plotting of calibration curves detract from the pedagogical value of
39 spectrophotometry. The project objective was a deliberately limited functionality that focused on
40 students learning the experimental aspects of spectrophotometry.
41
42
43
44
45
46
47

48 Other manufacturers offer small format spectrophotometry equipment that can address the space
49 considerations. For example, the Pasco wireless spectrometer (model PS-2600) is very small.
50 However, it requires a connected computer for the software that records the data or the company's
51 SPARK LXi Datalogger (model PS-3600). Likewise, Vernier has several compact options, the most
52 appropriate of which is the Go Direct® SpectroVis® Plus Spectrophotometer. The spectrophotometer
53
54
55
56
57
58

1
2
3 can be connected to a computer for recording data or to a Vernier LabQuest2 data-collection system.
4
5 Approaching this project with a Pasco or Vernier spectrophotometer connected to a computer doesn't
6
7 change the space issues significantly and would have required purchasing the spectrophotometers.
8
9 Using the SPARK LXi or LabQuest2 system would have been extremely space-efficient. Ultimately this
10
11 60 option was not selected, largely for economic considerations; the complete spectrometer and data-
12
13 collection system would be several times more expensive than the chosen solution.
14

15 **EQUIPMENT**

16 The spectrograph used for recording visible spectra in this work was a StellarNet Black Comet
17
18 C-25 (StellarNet Inc, Tampa, Florida). The particular model had a 2048-pixel CCD, 16-bit ADC, and a
19
20 65 $\Delta\lambda_{\text{effective}}$ of about 1 nm across the spectral range 179 to 1016 nm. The project was also tested on an
21
22 Ocean Optics HR4000 (Ocean Insights Inc., Largo, Florida). The software interface was developed in
23
24 Python (version 2.7 and 3.7, Python Software Foundation) on a Pi-top CEED (Pi-top, London, UK). A
25
26 wireless compact keyboard (model RT-MWK12+, Riitek Technology, Shenzhen, China) was selected to
27
28 minimize space requirements.
29

30 70 There are two versions of the project. One simply controls the spectrograph and would require
31
32 the use of existing light sources and cuvette holders. The second integrates a light source into the
33
34 software and Raspberry Pi hardware. Two light sources were investigated: a "cool white" LED (LTW-
35
36 2R3D7, Lite-On Inc.) and a "warm white" LED (334-15/X1C5-1QSA, Everlight Electronics Co Ltd).
37
38 The spectra of the two LEDs are shown in Figure 1 under conditions that place the maximum emission
39
40 75 at the upper limit of the spectrometer range. The choice between the LEDs is up to the user
41
42 depending on the spectral range in which they are interested. A cuvette holder was manufactured by
43
44 the author from Delrin thermoplastic for direct attachment to the face of the spectrograph. The
45
46 cuvette holder is designed for use with the integrated LED light source; other cuvette holders can be
47
48 used if an external light source is used. Fabrication plans for the cuvette holder are available in the
49
50 80 supplementary material. In the project version with an integrated LED light source, the brightness of
51
52 the LED is adjustable in the software interface via a pulse width modulation (PWM) circuit shown in
53
54 Figure 2. The PWM frequency can be adjusted in the code but is set initially to 2000 Hz which is
55
56 sufficiently fast to avoid the appearance of signal flickering at integration times greater than or equal
57
58

to 3 ms (the spectrometer minimum integration time). The PWM duty cycle was tested from 1 % to 100 % and found to produce constant LED intensity at all values.

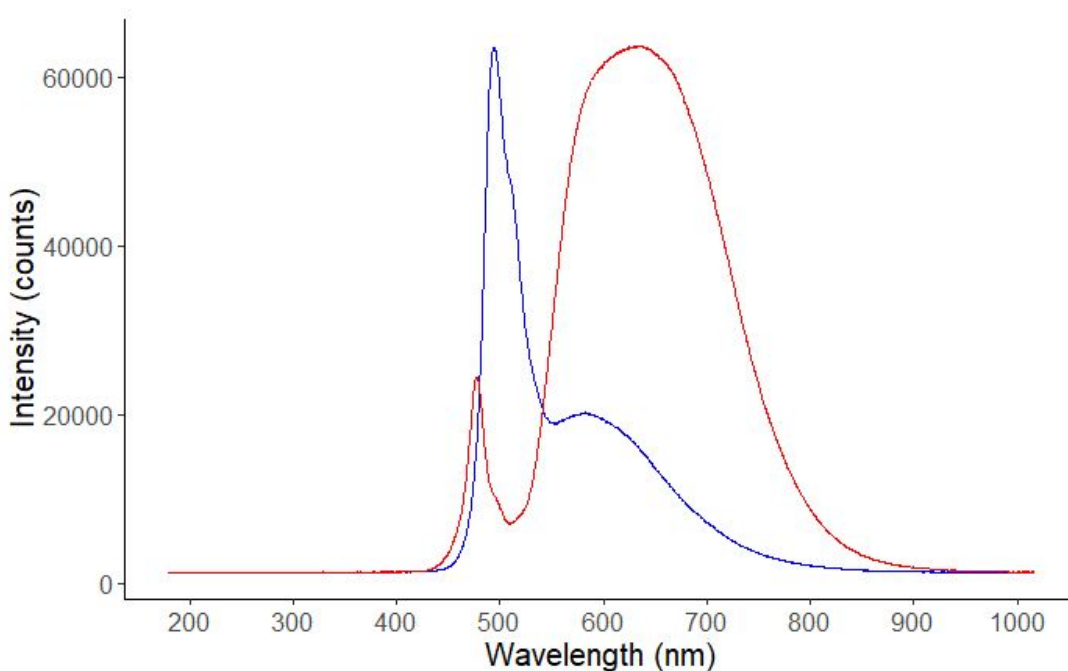


Figure 1. The emission spectra of the LEDs. The spectrum of the LTW-2R3D7 is drawn in blue and the 334-15/X1C5-1QSA is drawn in red. The spectra were recorded with 5-ms integration time and an LED duty-cycle of 99 % (red trace) and 31 % (blue trace) so that the maximum intensity values are equivalent and near the maximum of the spectrometer range.

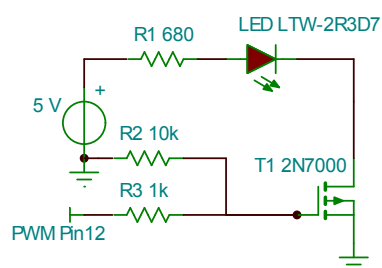


Figure 2. The pulse width modulation circuit for controlling the brightness of an LED light source. The pin number is relevant to a Raspberry Pi implementation of the circuit.

HAZARDS AND SAFETY PRECAUTIONS

There are not any inherent chemical hazards in the implementation of this project as a spectrophotometer. Its use in a teaching laboratory will need to respect the risks and hazards of the chemical system that is studied with the spectrophotometer.

RESULTS

The major product of the project is a simple software interface for performing spectrophotometry measurements. The software written for this project consists of a Python module that collects the physical parameters of the spectrograph at startup and provides an interactive interface for spectrophotometry. The code for the interface is located on GitHub (<https://github.com/acpo/>) with versions for StellarNet and Ocean Optics spectrographs. Instructions for installing the project are available in the supplemental information. There are versions available with and without an integrated LED light source. The integration of the LED light source provides control of all parts of the spectrophotometry instrument within the interface. A description of the user interaction with the interface is shown in Table 1 for the cases in which most experimental conditions are preset in the code (*e.g.*, for the least experienced users) and the case in which the user must work with all of the experimental conditions.

Table 1. User Interaction with for a Typical Spectrophotometric measurement

	Case: some conditions preset in code	Case: all conditions set by user
Select integration time such that 100 % T signal is on scale ^a (<i>text box</i>)	Not required	Required
Select wavelength boundaries appropriate to the light source (<i>text box</i>)	Not required	Required
Select the number of spectra to average ^b (<i>text box</i>)	Not required	Required
Record 0 % T (<i>button click</i>)	Required	Required
Record 100 % T (<i>button click</i>)	Required	Required
Switch to Absorbance mode (<i>button click</i>)	Required	Required
Select wavelength to monitor (<i>text box</i>)	Required	Required
Observe absorbances	Required	Required

^a if the integrated LED is used then the LED duty-cycle and the integration time both must be adjusted to give an on-scale signal.

^b averaging performs the arithmetic mean of the number of successive spectra indicated by this value, and acquisition requires a total time of the integration time × averages

An example of the complete system hardware and software interface are shown in Figure 3.

The Raspberry Pi-based system used in this project was useful to demonstrate that a very low-power computer was suitable. On the Pi-top system the response times of the interface and spectrum redrawing were fast enough for most students to be comfortable using it. The project could readily be

adapted to a small (3.5 or 5-inch) touch screen display to radically reduce the size of the device,
including the elimination of the keyboard.

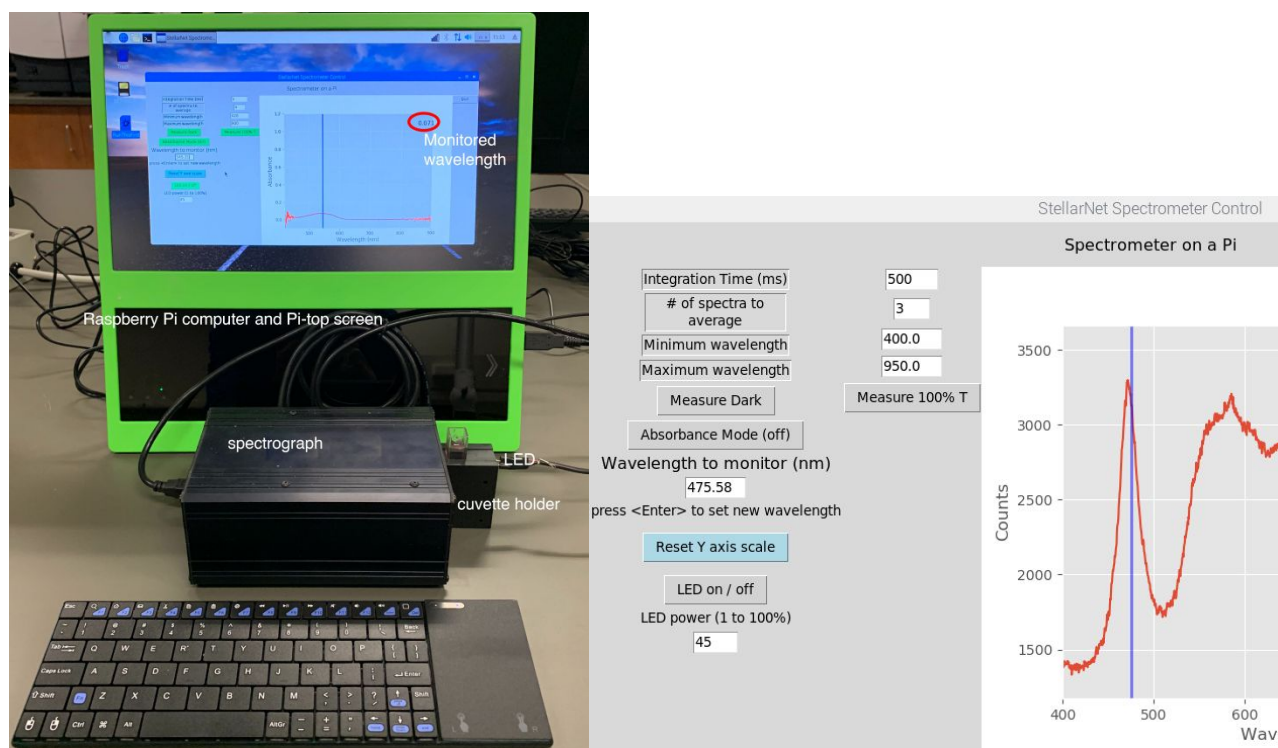


Figure 3. The complete spectrophotometer implementation (left) and the Python interface (right).

The project as implemented on a Pi-top computer and StellarNet spectrograph was tested with a quantitative analysis laboratory of 45 students divided into sections of about 12 each. Each student worked independently on a spectrophotometric determination of iron with *o*-phenanthroline by external calibration.^{10, 11} Relative to the student results of prior years in which a Spectronic 20 was used in this experiment, the data showed similar linearity in the range of 0 to 0.6 absorbance units. Linearity of data above 0.6 absorbance was slightly better with the compact spectrograph. It should be noted that the spectrophotometric performance of the system is largely dependent upon the particular hardware that a user selects, so the author's experience may not be the same as other users. However, in the transition from a Spectronic 20 ($\Delta\lambda_{\text{effective}} \approx 20$ nm) to a modern compact spectrograph ($\Delta\lambda_{\text{effective}}$ typically less than 2 nm) it would be unusual to find poorer performance related to deviations from the Beer-Lambert law.

The student experience with this project was similar to the Spectronic 20, with some advantages. Both systems required the student to record a 0 % and 100 % transmittance

1
2
3 measurement at the start of analysis. The Spectronic 20 accomplished this task via two analog dials
4
5 whereas in this project students clicked two buttons on the interface; effectively the student
6
7 experience is identical. Once the prerequisites were recorded, the absorbance at the selected monitor
8
9 wavelength was continuously displayed on the spectrum window (see Fig. 3). Selection of the
10
11 135 monitored wavelength was done via a numerical entry box on the interface. An important advantage
12
13 in this project was that the absorbance spectrum was displayed continuously during use. The
14
15 continuous display of the complete absorbance spectrum offers the instructor an opportunity to have
16
17 students easily explore the concept of λ_{\max} and other properties of the Beer-Lambert law.

18
19 A pedagogical advantage to this project is that the Python code for this project is available to be
20
21 140 edited and customized to an instructor's needs. The distributed version sets the wavelength range to
22
23 the range reported by the spectrograph firmware. An instructor has the choice to modify two lines of
24
25 code to preset the boundaries or to ask students to select with the interface the appropriate
26
27 boundaries for the light source that they are using. For example, if a white light LED is used, the
28
29 range can be restricted to 450 to 850 nm. Likewise, the integration time of the detector can be preset
30
31 145 in the code by the instructor to appropriately utilize the dynamic range of the detector or students can
32
33 be asked to explore the how the integration time can be used to make best use of the available light.
34
35 Additionally, the interface developed in this project gives the user access to the raw signal from the
36
37 spectrograph and can be used for emission (*e.g.*, fluorescence or chemiluminescence) measurements
38
39 which were never part of the Spectronic 20 scope of functions.

40
41 150 Lastly, to the author's knowledge, the interface developed in this project is the first free
42
43 spectrophotometry interface for StellarNet spectrometers on the Linux platform. While this project
44
45 was being tested, the company introduced a Mac OS compatible software, but there are not currently
46
47 any solutions for the Raspberry Pi or other Linux distributions. This project also offers a free, open-
48
49 source Linux interface to Ocean Optics instruments which the company does not currently provide.

51 155 **ASSOCIATED CONTENT**

52 Supporting Information

53 The Supporting Information is available on the ACS Publications website at DOI:

54
55 10.1021/acs.jchemed.XXXXXXX. [ACS will fill this in.]
56
57
58

Notes to guide the installation of the project. NotesForInstallation.pdf (PDF)

Cuvette holder fabrication plans. FabricationPlans.pdf (PDF)

A description of each functional unit in the software. SoftwareDescription.pdf (PDF)

A description of the hardware requirements for Raspberry PI, Windows, and Apple systems.

HardwareDescription.pdf (PDF)

AUTHOR INFORMATION

Edward C. Navarre

*E-mail: enavarr@siue.edu

ACKNOWLEDGMENTS

The students of the Summer and Fall 2019 sessions of the quantitative analytical chemistry laboratory at SIUE are thanked for their patience in testing this project. The reviewers of the manuscript are thanked for their comments that helped to reveal several issues with earlier versions of the project.

Edward C. Navarre received funding from the Southern Illinois University Edwardsville Excellence in Undergraduate Education Program for part of this project.

REFERENCES

1. Brown, T. L. *Chemistry: The Central Science*, 13th edition.; Pearson: Boston, 2015. pp. 1019-1020.
2. Robinson, J. K.; McMurry, J.; Fay, R. C. *Chemistry*, 8th Ed.; Pearson Education, Inc: Hoboken, NJ, 2020. pp. 929-930.
3. Skoog, D. A. *Fundamentals of Analytical Chemistry*, 9th Ed.; Cengage - Brooks/Cole: Belmont, CA, 2012. pp. 650-674.
4. Harris, D. C.; Lucy, C. A. *Quantitative Chemical Analysis*, 9th Ed.; W.H. Freeman & Company: New York, 2016. pp. 432-443.
5. Griffiths, P. R. Whither "Quant"? An Examination of the Curriculum and Testing Methods for Quantitative Analysis Courses Taught in Universities and Colleges in the Western USA. *Anal. Bioanal. Chem.* **2008**, 391 (3), 875–880. <https://doi.org/10.1007/s00216-008-2080-2>.
6. Grasse, E. K.; Torcasio, M. H.; Smith, A. W. Teaching UV–Vis Spectroscopy with a 3D-Printable Smartphone Spectrophotometer. *J. Chem. Educ.* **2016**, 93 (1), 146–151. <https://doi.org/10.1021/acs.jchemed.5b00654>. And references 3 – 9 therein.

-
- 1
2
3 7. Vanderveen, J. R.; Martin, B.; Ooms, K. J. Developing Tools for Undergraduate Spectroscopy: An
4 Inexpensive Visible Light Spectrometer. *J. Chem. Educ.* **2013**, *90* (7), 894–899.
5 190 <https://doi.org/10.1021/ed300396x>.
6
7 8. Bougot-Robin, K.; Paget, J.; Atkins, S. C.; Edel, J. B. Optimization and Design of an Absorbance
8 Spectrometer Controlled Using a Raspberry Pi To Improve Analytical Skills. *J. Chem.*
9 *Educ.* **2016**, *93* (7), 1232–1240. <https://doi.org/10.1021/acs.jchemed.5b01006>.
10
11 9. Cell Phone Spectrometer
12 195 [https://www.asdlib.org/onlineArticles/elabware/Scheeline_Kelly_Spectrophotometer/](https://www.asdlib.org/onlineArticles/elabware/Scheeline_Kelly_Spectrophotometer/index.html)
13 [index.html](https://www.asdlib.org/onlineArticles/elabware/Scheeline_Kelly_Spectrophotometer/index.html) (accessed Oct 28, 2019).
14
15 10. Fortune, W. B.; Mellon, M. G. Determination of Iron with O-Phenanthroline: A
16 Spectrophotometric Study. *Ind. Eng. Chem. Anal. Ed.* **1938**, *10* (2), 60–64.
17 <https://doi.org/10.1021/ac50118a004>.
18 19 200
19 20 11. Atkins, R. C. Colorimetric Determination of Iron in Vitamin Supplement Tablets. A General
21 Chemistry Experiment. *J. Chem. Educ.* **1975**, *52* (8), 550. <https://doi.org/10.1021/ed052p550>.
22
23
24
25
26
27
28
29
30
31
32
33
34
35
36
37
38
39
40
41
42
43
44
45
46
47
48
49
50
51
52
53
54
55
56
57
58